Instruction Manual

MEGLY

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Disclaimer

Our Company is not responsible for damage resulting from fire, earthquakes, riots, terrorism or other events; from intentional or unintentional misuse by the customer; or from abnormal conditions.

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ENGLISH

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
Safety Precautions

- This product has built-in magnets in the locations indicated in the figure.

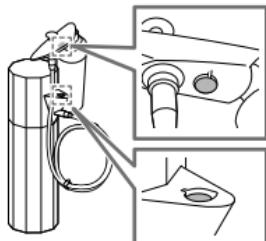Never use this product together with the following medical electrical apparatuses. Accidents or medical electrical apparatus malfunction may result, causing serious physical disabilities.

- Implantable medical electrical apparatuses such as pacemakers
 - Variable pressure cerebrospinal fluid shunts*
- * A medical apparatus embedded under the scalp that is used for a long time period to bypass cerebrospinal fluid into the abdominal cavity and is capable of changing the flow rate of cerebrospinal fluid by changing the valve pressure magnetically.

- Medical electrical apparatuses for life support such as a heart-lung machine
- Wearable medical electrical apparatuses such as an electrocardiograph
- Other medical electrical apparatuses that are affected by magnetic forces
- Remove items you are wearing such as jewelry, metal hairpins, necklaces, pierced earrings, glasses, hearing aids, earphones, or headphones. They may be attracted by the magnets of this product, causing damage or unexpected injuries.
- Never use this product near metal products. This product may be attracted, causing damage or unexpected injuries.
- Do not place clocks, magnetic cards, personal computers, electronic clocks, magnetic disks, cell phones, smartphones, tablets, music players, remote controllers, or other items that are affected by magnetic forces near the magnets of this product. There is a risk of damage to equipment or data.
- The following persons should not use this

product. There is a risk of an accident or problems to the body or skin.

- Persons prone to metal allergies
- Persons with sensory abnormalities of the skin
- Persons with infected or inflamed skin
- Small children
- Persons with physical disabilities
- Persons who have been drinking alcohol
- Persons in poor health
- Anyone with the following conditions should not use the product without first consulting a physician. There is a risk of an accident or problems to the body or skin.
 - Persons prone to allergies
 - Persons with sensitive skin
 - Persons with atopic dermatitis
 - Persons who have experienced skin rash with cosmetics or disinfectants
 - Persons with malignant tumors
 - Persons with acute illnesses
 - Persons with severe blood circulation disorders or blood diseases

- Persons who are pregnant, may be pregnant, or are in their menstrual cycle (offset hormone balance can cause skin problems.)

- Persons receiving treatment from a medical institution

* Even if none of the conditions above apply to you, please consult with your physician if you have any concerns about using the product.

▪ Do not use the product on the following body parts. There is a risk of an accident or problems to the body or skin.

- On wounds, or mucous membranes such as the eyes or mouth
- On skin abnormalities (swelling, pimples, a festering condition, or severe sunburn)
- On body parts where hair removal treatment was performed the same day
- On body parts where plastic surgery has been performed
- On body parts with metal, plastic, or silicone implants

* Please do not use anywhere other than as specified in this manual.

- 4
- Due to the characteristics of CO₂ gas, the area of the skin on which the gas is sprayed may redden for some people. If you have any concerns about using this product on your skin, try it out on an area of the body that is less noticeable to make sure that your skin does not react adversely.
 - When using the product, make sure your skin is free of irregularities. Immediately discontinue use if a cosmetic causes skin irritation, more specifically, if any of the following occurs. If you continue using the cosmetic, the symptoms may worsen, so consult a dermatologist.
 - (1) Redness, swelling, itching, irritation, color loss (such as white spots), or darkening
 - (2) Conditions such as those described above that develop when the skin is exposed to direct sunlight after using the product
 - Avoid inhaling a lot of CO₂ gas and cosmetics when spraying around your nose or mouth.
 - Do not spray directly into eyes, nor rub eyes right after use.
 - Close your eyes when spraying around the area of your eyes. If you get any spray in your eyes, rinse your eyes immediately with water. If you experience any unusual symptoms, seek medical assistance.
 - Do not use in any way except for its intended purpose. There is a risk of an accident, injury, damage, or gas leak.
 - Do not disassemble or modify this product. There is a risk of an accident, injury, or damage.
 - Use the Tube, Handpiece, and other parts that are exclusively for this product. Using other parts may cause accidents, problems on your skin, injury, or damage.
 - Do not pull the Handpiece hard while using. The main unit may fall and result in an accident, injury, or damage.
 - Do not disconnect the Tube from the tube joint until the CO₂ in the CO₂ gas cartridge is fully consumed (when using or storing). There is a risk of an accident, injury, damage, or gas leak.
 - Make sure that the Tube is free of crimping, crushing, or damage during use or storage. There is a risk of an accident, injury, damage, or gas leak.
 - Never use or store this product in the following locations. There is a risk of an accident, injury, breakage, or damage.
 - Places prone to vibration or impacts and unstable locations

- Places within the reach of children (there is a risk of accidental swallowing of small parts such as the Cleaning Wire and the Tube Clips)
 - Places exposed to direct sunlight or environmental temperatures of less than 5°C (41°F) or over 40°C (104°F), such as near an open flame or in a vehicle
 - Places exposed to high temperatures such as on electric carpets, or near heating appliances or lighting devices that produce heat
 - Places where there is a large atmospheric pressure or temperature change such as at high altitudes
 - Places exposed to water, or places with a lot of humidity such as a bathroom
 - Places prone to corrosion by a corrosive substance such as salt or a chemical agent
 - Dusty areas
 - Enclosed spaces
- Do not move this product when cosmetics are adhered to it or your hands. The product may slip and drop.
 - When you move while holding this product,

- hold the Main Unit (Lower) and keep it upright, taking care not to drop it. Dropping it may result in an accident or injury as well as damage to this product.
- Take care not to prick fingers or any other body parts with the tip of the Cleaning Wire.

Handling of CO₂ Gas Cartridge

- Only use MEGLY specified CO₂ gas cartridges marked "Carbonated Mist Unit Use Only". Using other cartridges or parts may cause accidents, problems on your skin, injury, or breakage.
- The CO₂ gas cartridge is filled with high pressure gas. Thus you cannot carry it onto an airplane in a carry-on nor checked baggage, according to the Civil Aeronautics Law. In addition, you cannot send it by airfreight.
- Do not use a CO₂ gas cartridge that shows any signs of rust, scratching, or deformation. There is a risk of rupture or explosion.
- Do not drop the CO₂ gas cartridge or subject it to shock. There is a risk of an accident, injury, or damage.
- Do not discard or leave a CO₂ gas cartridge that is not empty lying around. There is a risk of rupture or explosion.

- 6
- Do not cut or bore holes in a CO₂ gas cartridge with a method other than the correct attaching method described in P.12 "Attach the CO₂ gas cartridge". There is a risk of rupture or explosion.
 - Before attaching a CO₂ gas cartridge, wipe any foreign matter (such as hair and dust) from the tip using a dry cloth or other material.
 - The CO₂ gas cartridge may become extremely cold when used continuously for 20 seconds or more, or when repeatedly used in a short time. In addition, if it is not attached to the Main Unit correctly and CO₂ gas is continuously sprayed, the CO₂ gas cartridge may also become extremely cold. Follow the precautions below when using the product. There is a risk of damage, gas leak, or injury (such as frostbite).
 - * If you accidentally touch the cold CO₂ gas cartridge and have skin irritation, do not rub the affected region, but gently apply lukewarm water, lightly wrap the area in gauze, and immediately seek medical assistance.
 - Do not use the product continuously for more than 20 seconds.
 - Do not repeatedly use the product in a short time. After 1 use, wait 5 minutes or longer before the next use.
 - When attaching or removing the CO₂ gas cartridge, do not touch it with bare hands or wet hands. Always make sure to wear nonslip gloves or rubber gloves.
 - When screwing the cartridge into the attachment socket of the Main Unit (Upper), take care not to hit the Main Unit (Upper) or its surroundings. There is a risk of an accident, injury, or damage.
 - Attach the CO₂ gas cartridge to the Main Unit (Upper) in a straight, vertical position. Do not screw it in diagonally.
 - Do not loosen the cartridge while screwing it in. There is a risk of a gas leak.
 - Ignore minor gas discharge with a "hissing" sound while attaching, and screw in the cartridge in a single quick motion. If you suspend it, gas may discharge and result in injury or damage.
 - After you feel resistance, screw the cartridge with force in a single quick motion until it no longer turns.

- Make sure that there is no leaking gas sound after screwing it in.
- When you hear a CO₂ gas leak sound after attachment, retighten the cartridge. If it cannot be retightened, leave it outside until the gas leak sound stops, and remove it.
- Do not contact with the gas being sprayed or the CO₂ gas cartridge during use. Do not direct the gas at a person.
- If the tube joint is removed while CO₂ gas remains, the gas is discharged with great force. If it occurs, while taking care not to touch the CO₂ gas cartridge with bare hands, leave it outdoors until the gas is fully discharged.
- After attaching the CO₂ gas cartridge, do not remove it until all the CO₂ gas has been used.
- After attaching the CO₂ gas cartridge, make sure to always use and store this product in an upright position. Using it sideways may cause liquefied CO₂ gas to flow into the Main Unit (Upper) which could result in damage to the product, a gas leak, or injury.
- Never use or store CO₂ gas cartridges in the following locations. There is a risk of an accident, injury, breakage, or damage.
 - Places prone to vibration or impacts and unstable locations
 - Places within the reach of children
 - Places exposed to direct sunlight or environmental temperatures of less than 5°C (41°F) or over 40°C (104°F), such as near an open flame or in a vehicle
 - Places exposed to high temperatures such as on electric carpets or near heating appliances or lighting devices that produces heat
 - Places where there is a large atmospheric pressure or temperature change such as at high altitudes
 - Places exposed to water, or places with a lot of humidity such as a bathroom
 - Places prone to corrosion by a corrosive substance such as salt or a chemical agent
 - Dusty areas
 - Enclosed spaces

Precautions for Use

- Place the Handpiece on the Main Unit and store the Main Unit in an upright position when suspending or after use. If the exclusive cosmetic bottle falls, liquid may leak or this product may be damaged.
- The spray amount of CO₂ gas or cosmetics may change depending on the environmental temperature. This is due to the structure of this product and not a problem with the product.
- Clogging may occur depending on the usage environment and frequency, and may prevent cosmetics from being sprayed. See P.19 "When removing clogging of the cosmetic" to clean the product periodically.
- If a cosmetic adheres to the back side of the Exclusive Cosmetic Cap, it may cause spraying failure due to clogging (e.g. spraying CO₂ gas without cosmetics). To prevent cosmetic clogging, follow the cautions below.
 - Do not lay the bottle down; do not tilt it more than necessary
 - Do not shake the bottle violently
- Do not use the bottle and the Handpiece while turning them upside down
- Do not drop the bottle
- Remove the exclusive cosmetic bottle and put the cap supplied with the exclusive cosmetic bottle back on to store it when the product remains unused for a long period (approximately 1 week or more). See P.24 "For long-term storage" when the product remains unused for a prolonged period (approximately 1 month or more).

Parts Names, Accessories

9

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
ENGLISH

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

Main Unit (Upper)

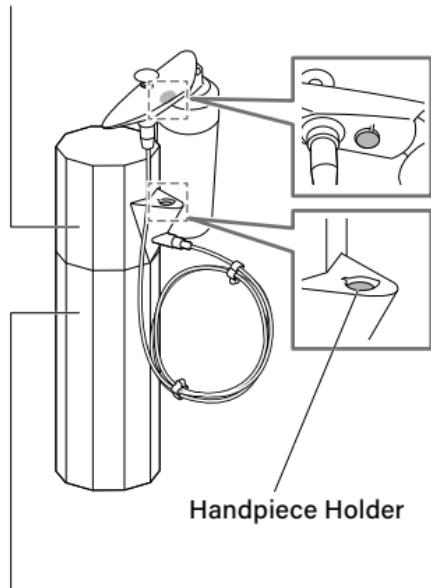
Main Unit (Lower)

Spray Button

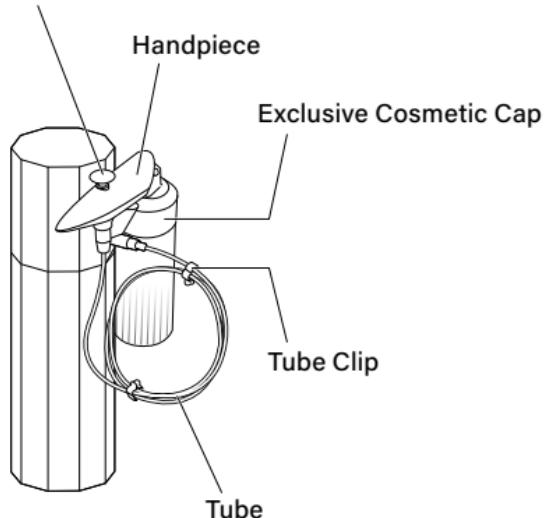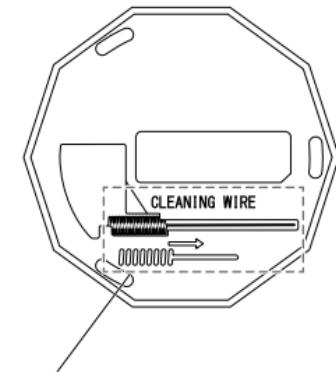Bottom

Cleaning Wire Housing
(P.24 "Storing the Cleaning Wire")

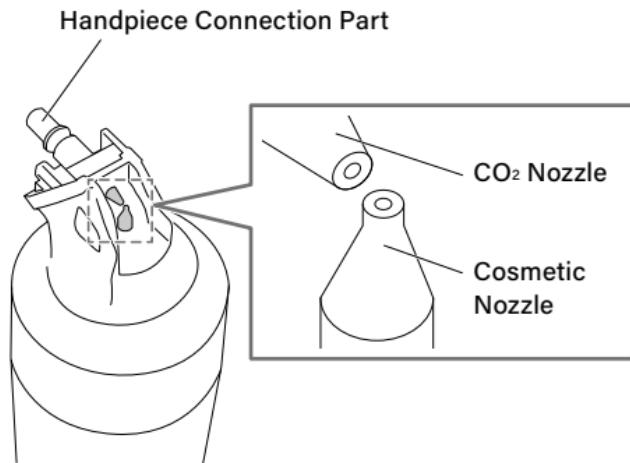**Accessories**

Cleaning Wire (stored inside the Cleaning Wire Housing at the time of purchase)	1
Instruction Manual (this document)	1
Warranty Certificate (printed on the Instruction Manual case)	1

See the MEGLY brand site for details.

Exclusive CO₂ Gas Cartridge

Model Number: GS-AB-00A (1)

PR-GC2, E1025GC-02 (2)

PR-GC6, E1025GC-06 (6)

Exclusive Cosmetics

This product supports only the "MEGLY" brand cosmetics. Using other brand cosmetics may result in clogging, spraying failure, or damage.

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ENGLISH

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
Preparation

Attach the CO₂ gas cartridge


- Read P.5 "Handling of CO₂ Gas Cartridge" before attaching the CO₂ gas cartridge.
- A CO₂ gas cartridge for the first use comes with the kit for the customers who purchase the MEGLY start kit. See P.11 "Optional Accessories" to purchase more for subsequent use.

- Only use MEGLY specified CO₂ gas cartridges marked "Carbonated Mist Unit Use Only". See P.11 "Optional Accessories" for details.
- If it is not attached to the Main Unit correctly and CO₂ gas is continuously sprayed, the CO₂ gas cartridge may become extremely cold. Do not touch it with bare hands, but always wear nonslip gloves or rubber gloves when handling it.


1. Holding the Main Unit (Lower), turn the Main Unit (Upper) in the direction of the arrow.

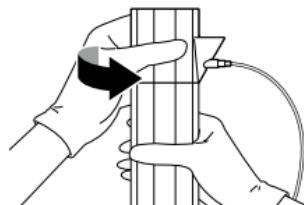2. Remove the Main Unit (Upper) in the direction of the arrow.

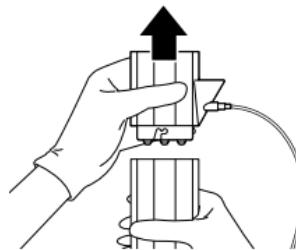
3. Hold the CO₂ gas cartridge vertically and press it against the attachment socket.

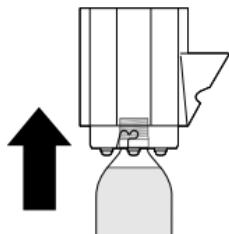- Before attaching a CO₂ gas cartridge, wipe any foreign matter (such as hair and dust) from the tip using a dry cloth or other material.
- Do not attach the CO₂ gas cartridge in a downward direction or sideways.
- Attach the CO₂ gas cartridge to the Main Unit (Upper) in a straight, vertical position. Do not screw it in diagonally.

4. Screw in the CO₂ gas cartridge in the direction of the arrow.



- After you feel resistance, screw the cartridge with force in a single quick motion until it no longer turns.

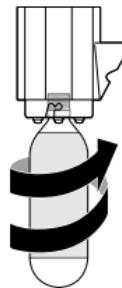- Do not loosen the cartridge while screwing it in. There is a risk of a gas leak.
- Ignore minor gas discharge with a "hissing" sound while attaching, and screw in the cartridge in a single quick motion.
- Make sure that there is no leaking gas sound after screwing it in.

- 14
5. Align the projection part of the Main Unit (Upper) and the recessed part of the Main Unit (Lower) and insert the projection into the recess.
 6. Turn the Main Unit (Upper) in the direction of the arrow.

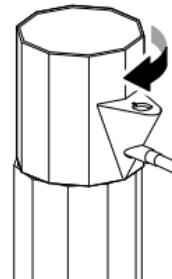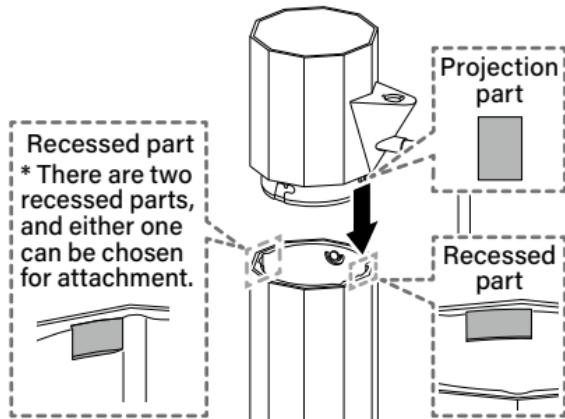Attach an exclusive cosmetic


1. Remove the cap of the exclusive cosmetic.

- Do not discard the cap because it is necessary for storing the cosmetic.
- An exclusive cosmetic for the first use comes with the kit for the customers who purchase the MEGLY start kit. See P.11 "Optional Accessories" to purchase more for the subsequent use.

2. Attach the Exclusive Cosmetic Cap for this product to the exclusive cosmetic.

3. Insert the connection part of the Exclusive Cosmetic Cap into the mounting hole on the Handpiece.

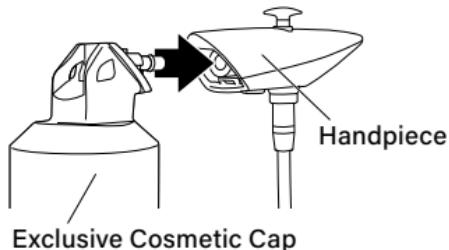- Insert it fully until it clicks. If it is loose, the cosmetic may not be sprayed or the cap may come off during use.


How to Use

Check the spraying of mist

After lightly shaking the exclusive cosmetic bottle sideways, press the Spray Button to check the spraying of mist.

- Aim the CO₂ Nozzle of the Handpiece at the back of your hand and press the Spray Button to confirm that mist is sprayed.
- Make sure that the Tube is free of crimping or crushing.

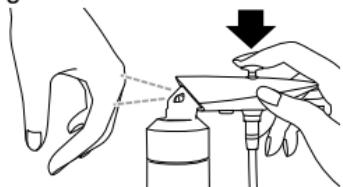* If no cosmetic is sprayed, the Tube may be clogged with the cosmetic. Follow the steps in P.19 "When removing clogging of the cosmetic" to clean the Tube. If cleaning cannot solve the spraying failure, or if clogging occurs frequently, please contact MEGLY Customer Support.

Spray the cosmetic on the face

Taking care so that the cosmetic does not get into your eyes, spray an appropriate amount on your entire face evenly.


- Always keep the unit 5 to 10 cm away from the area you are spraying.
- Avoid inhaling a lot of CO₂ gas and cosmetics when spraying around your nose or mouth.
- When using this product for the first time, spray the cosmetic on an area of the body that is less noticeable to make sure that your skin does not react adversely.
- To prevent falling, place the Main Unit on a stable place or hold it in your hand during use.
- Do not use the product with the Main Unit in the upside-down or sideways position.

Replace the CO₂ gas cartridge


- If CO₂ gas is not released when the Spray Button is pressed, the gas cartridge is empty. In this case, replace the empty cartridge with a new one by following the steps below.
- Read P.5 "Handling of CO₂ Gas Cartridge" before replacing the CO₂ gas cartridge.

- After attaching the CO₂ gas cartridge, do not remove it until all the CO₂ gas has been used.
- The CO₂ gas cartridge may become extremely cold if it is removed immediately after use or while CO₂ gas remains. Do not touch it with bare hands, but always wear nonslip gloves or rubber gloves when handling it.


- Press the Spray Button to confirm that the gas has been completely exhausted.

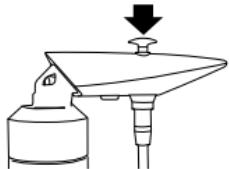
- Hold the Main Unit (Lower) tightly, and turn the Main Unit (Upper) in the direction of the arrow.

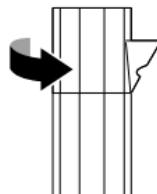- Pull out the Main Unit (Upper) in the direction of the arrow, and remove it with the CO₂ gas cartridge.

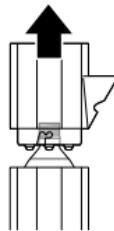4. Hold the Main Unit (Upper) tightly, and turn the CO₂ gas cartridge in the direction of the arrow and remove it.

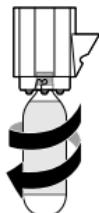- Do not remove the CO₂ gas cartridge in a downward direction or sideways.

5. Attach the replacement CO₂ gas cartridge.

- See P.12 "Attach the CO₂ gas cartridge".

When removing smudges attached to the surface

- Do not wipe the product with paint thinner, benzine or alcohol. There is a risk of product damage such as deformation or cracking of parts.


1. Wipe off smudges on the Main Unit with a dry soft cloth.

- If smudges do not come off, wipe them off with a damp, tightly wrung out cloth.

When removing clogging of the cosmetic


If only CO₂ gas is sprayed but no cosmetic is sprayed, clean the product by following steps ①, ②, and ③ below.

- Cleaning the Cosmetic Nozzle(P.20)
- Cleaning the CO₂ Nozzle(P.21)
- Soaking the Exclusive Cosmetic Cap for cleaning(P.21)

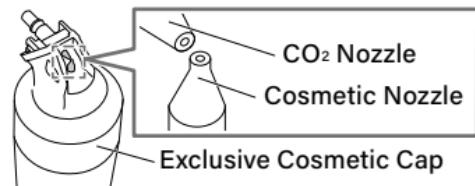Use the Cleaning Wire for cleaning.

Bottom of the Main Unit

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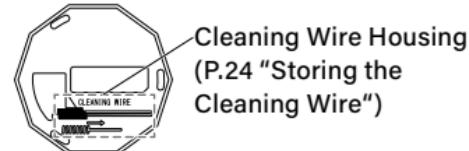ENGLISH

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
① Cleaning the Cosmetic Nozzle

1. Remove the Handpiece from the Exclusive Cosmetic Cap.

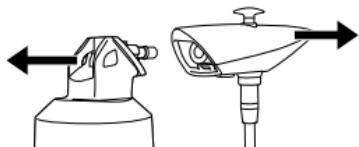2. Remove the exclusive cosmetic bottle from the Exclusive Cosmetic Cap.

- Put the cap supplied with the exclusive cosmetic bottle back on to store the removed exclusive cosmetic.

3. Insert the included Cleaning Wire straight from the center hole on the back of the Exclusive Cosmetic Cap, and move it up and down 2 or 3 times to push out the clog.

- Move the Cleaning Wire so that its tip comes out of the Cosmetic Nozzle on the opposite side of the Exclusive Cosmetic Cap.

- The Cleaning Wire is inside the Cleaning Wire Housing at the time of purchase. See P.24 "Storing the Cleaning Wire" for storage and extraction.

- Take care not to prick fingers or any other body parts with the tip of the Cleaning Wire.
- Handle the Cleaning Wire with care as it is thin and easily bent.

* If the steps of "① Cleaning the Cosmetic Nozzle" cannot solve the spraying failure, try "② Cleaning the CO₂ Nozzle".

② Cleaning the CO₂ Nozzle

1. Remove the Handpiece and the exclusive cosmetic bottle from the Exclusive Cosmetic Cap.
- This procedure is unnecessary when proceeding to the CO₂ Nozzle after "① Cleaning the Cosmetic Nozzle".
2. Insert the included Cleaning Wire straight from the hole on the Handpiece connection side, and move it up and down 2 or 3 times to push out the clog.


● Move the Cleaning Wire so that its tip comes out of the CO₂ Nozzle. The Cleaning Wire may touch the Cosmetic Nozzle depending on the angle. Do not press hard.

- * If the steps of "① Cleaning the Cosmetic Nozzle" and "② Cleaning the CO₂ Nozzle" cannot solve the spraying failure, try also "③ Soaking the Exclusive Cosmetic Cap for cleaning".


③ Soaking the Exclusive Cosmetic Cap for cleaning

1. Soak the Exclusive Cosmetic Cap in lukewarm water (approximately 35°C (95°F)) for approximately 5 minutes.

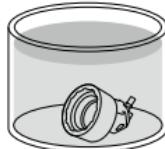2. Wash the back of the Exclusive Cosmetic Cap with tap water.


● Shake off excess water and dry in the shade.

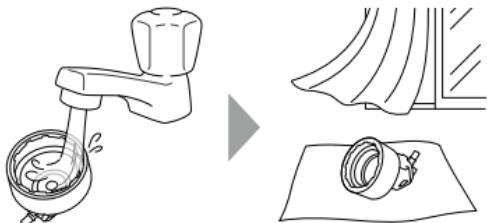* If these cleaning steps cannot solve the spraying failure, or if clogging occurs frequently, please contact MEGLY Customer Support.


22 When replacing cosmetics

When replacing exclusive cosmetic bottles,
wash the back of the Exclusive Cosmetic Cap
with tap water, shake off excess water, and dry
in the shade before replacement to prevent
mixing of cosmetics or clogging with cosmetics.

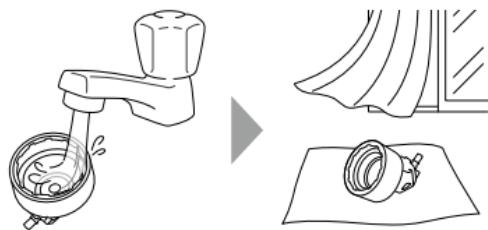
Storage of the Main Unit

- After placing the Handpiece in the position indicated in the figure, store the Main Unit on a stable surface, making sure that the Main Unit and the exclusive cosmetic bottle will not tip over.
- Tube length can be adjusted using the Tube Clips.

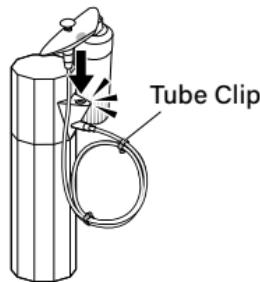- The Handpiece can be placed either on the right side or the left side.

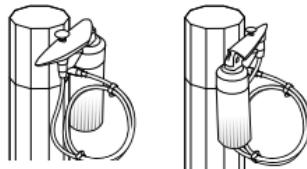- Never store or place CO₂ gas cartridges in the following locations. There is a risk of an accident, injury, breakage, or damage.
 - Places prone to vibration or impacts and unstable locations
 - Places within the reach of children (who may accidentally swallow small parts such as the Cleaning Wire and the Tube Clips)
 - Places exposed to direct sunlight or environmental temperatures of less than 5°C (41°F) or over 40°C (104°F), such as near an open flame or in a vehicle
 - Places exposed to high temperatures such as on electric carpets or near heating appliances or lighting devices that produces heat
 - Places where there is a large atmospheric pressure or temperature change such as at high altitudes
 - Places exposed to water, or places with a lot of humidity such as a bathroom
 - Places prone to corrosion by a corrosive substance such as salt or a chemical agent
 - Dusty areas
 - Enclosed spaces

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
ENGLISH

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- Make sure that the Tube is free of crimping, crushing, or damage during storage or placement. There is a risk of an accident, injury, damage, or gas leak.

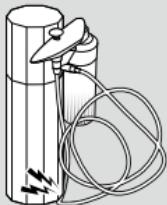For long-term storage

For long-term storage, make sure to follow the steps below before storage.

1. Remove the CO₂ gas.

- See P.27 "To discharge gas before disposal"

2. Remove the Exclusive Cosmetic Cap and clean the product.


- See P.19 "When removing clogging of the cosmetic"

- See the precautions provided on the cosmetic package or bottle for the storage method and storage life of the cosmetic.

Storing the Cleaning Wire


The Cleaning Wire can be stored on the bottom of the Main Unit (Lower).

Cleaning Wire Housing (Bottom of the Main Unit)

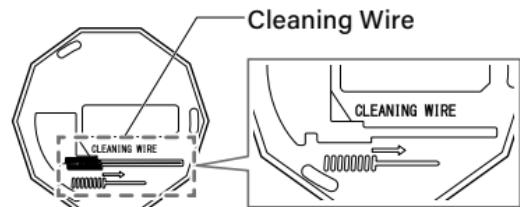
- The Cleaning Wire is inside the Cleaning Wire Housing at the time of purchase.
- Take care not to prick fingers or any other body parts with the tip of the Cleaning Wire.
- Handle the Cleaning Wire with care as it is thin and easily bent.

- * To prevent backflow or leakage of the cosmetic in the bottle when the Main Unit is turned upside down, please remove the exclusive cosmetic bottle from the Handpiece before storing the Cleaning Wire.

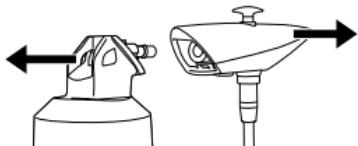
Storage

1. Insert the tip of the Cleaning Wire into (a) of the figure.
2. Hook the front of the spring part to (b), and push so that the back of the spring part is housed in (c).

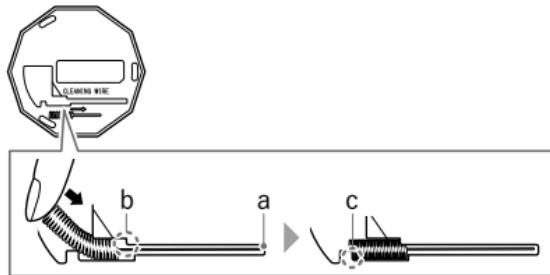Pulling out

1. Push the spring part from the back to compress it, and remove the Cleaning Wire from (c).
● If it is hard to push it with your fingers, use something such as a hairpin.
2. Pinch the back of the spring part, and pull it out.

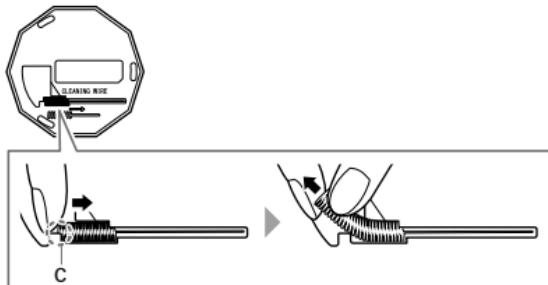
Disposal Method

Disposal of the Main Unit

- Dispose of the product according to the regulations of your local government.
- Discontinue use if cuts or cracks develop in this product, remove the CO₂ gas cartridge safely by referencing P.27 "To discharge gas before disposal", and dispose of the product.

Disposal of CO₂ gas cartridges

- Incinerating an unused CO₂ gas cartridge without a hole and filled with CO₂ gas can result in rupture or explosion.
- Used CO₂ gas cartridges have a hole on the cap (sealing cap). Check that the gas is used up by checking if there is no leak sound from the hole, and dispose of the cartridge according to the regulations of your local government (CO₂ gas cartridge body material: steel).

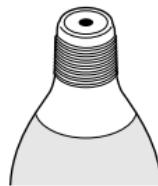
- To dispose of unused CO₂ gas cartridges, empty the gas by following the method below.

To discharge gas before disposal

- * Before attaching, read the precautions in P.5 "Handling of CO₂ Gas Cartridge".
- * Be sure to do it outside.
- * The CO₂ gas cartridge becomes extremely cold when the gas is discharged. Always wear nonslip gloves or rubber gloves when handling it.
- * Be sure to slowly and carefully loosen the CO₂ gas cartridge to prevent CO₂ gas from being discharged with great force.

1. Take this product outside before discharging the gas.
2. Attach the CO₂ gas cartridge from which you want to discharge gas to the Main Unit (Upper).
3. Slowly and carefully turn the CO₂ gas cartridge twice to loosen it, and check that there is a leak sound of the CO₂ gas.

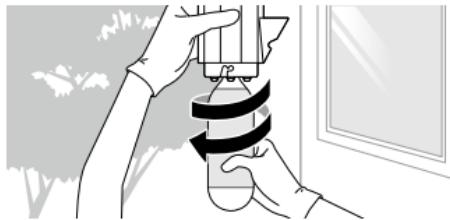4. After the leak sound stops, remove the CO₂ gas cartridge.

Troubleshooting

Read this section before sending out for repair.

Condition	Cause	Countermeasure
There is a hissing sound (CO ₂ gas leak sound).	The Tube, the Main Unit, and the Handpiece are not correctly connected.	Check that the Tube, the Main Unit, and the Handpiece are correctly connected. Retighten the tube joint if it is loose.
	The CO ₂ gas cartridge is not fully screwed in.	Fully screw in the CO ₂ gas cartridge to the Main Unit (Upper). See P.12 "Attach the CO ₂ gas cartridge".
CO ₂ gas is not released when the Spray Button on the Handpiece is pressed.	The CO ₂ gas cartridge is not fully screwed in.	Fully screw in the CO ₂ gas cartridge to the Main Unit (Upper). See P.12 "Attach the CO ₂ gas cartridge".
	The CO ₂ gas cartridge is empty.	Replace it with a new CO ₂ gas cartridge. See P.17 "Replace the CO ₂ gas cartridge".

Condition	Cause	Countermeasure
No cosmetic is sprayed.	The connection part of the Exclusive Cosmetic Cap is detached from the Handpiece.	Securely attach the connection part of the Exclusive Cosmetic Cap to the Handpiece.
	The cosmetic has caused clogging inside the Exclusive Cosmetic Cap.	Clean it by following the steps on P.19 "When removing clogging of the cosmetic".
The CO ₂ gas cartridge is empty even though it has not been used.	For safety reasons, a small amount of gas gradually leaks even if it is not used. The cartridge becomes empty because it is not used for a long period of time.	Replace it with a new CO ₂ gas cartridge. See P.17 "Replace the CO ₂ gas cartridge".

- * If the countermeasures fail to correct your problem or you encounter problems not covered in this section, immediately discontinue use of the product and contact MEGLY Customer Support.

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ENGLISH

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Q&A

Q1. Can any cosmetic lotions be used?

A1. This product supports only the cosmetics exclusively for MEGLY. Using other brand cosmetics may result in clogging, spraying failure, or damage.

Q2. How long can the CO₂ Gas Cartridge be used?

A2. The usable period is four to five weeks when it is used twice a day, for 10 seconds each time, although the actual spray amount depends on the usage environment. When gas is reduced earlier than usual, there is a possibility of a gas leak. Check if the CO₂ gas cartridge is fully screwed in.

Q3. What should I do if the cosmetic gets in my eyes?

A3. Rinse your eyes with water immediately. If you experience any unusual symptoms, seek medical assistance.

Product Specifications

31

Product name	MEGLY
Model number	GS-AA-02A
Product size	[Main Unit] Approx. 197 mm × 80 mm × 55 mm (excluding Handpiece) [Handpiece] Approx. 48 mm × 94 mm × 35 mm (including Spray Button and Exclusive Cosmetic Cap) [Tube] Approx. φ 4 mm × 1000 mm
Materials	[Main Unit] ABS, PC+ABS, aluminum, polyurethane, polyacetal, brass, polypropylene, PBT, stainless steel, neodymium magnet, steel, silicone rubber, nitrile rubber [Tube] Polyurethane, brass
Items included/ accessories	Items included: Main Unit (1) Accessories: Cleaning Wire (1), Instruction Manual (1), Warranty Certificate (printed on the Instruction Manual case) (1)
Weight	Approx. 315 g
Adjustment pressure	Approx. 0.15 MPa
Recommended usage temperature	5 to 33°C
Safety release	Functions on or below 0.9 MPa
Manufacturer	MTG Co., Ltd. (Japan) Country of origin: China

- Specifications may change without notice for further product improvements. Thank you for your understanding in advance.
- This product was planned, developed, and designed in Japan.

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ENGLISH

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Warranty and After-sales Service

Warranty certificate

Be sure to read the Warranty Certificate carefully and keep it in a safe place together with documents with the date of purchase clearly stated (such as the invoice).

Warranty period

The warranty period is 1 year from the date of purchase.

- * If there are any problems during use, contact MEGLY Customer Support.

Minimum retention period for service functional parts

The minimum retention period for service functional parts for this product is 5 years from the discontinuation of manufacture. "Service functional parts" means parts that are required to maintain the functionality of this product.

When requesting repairs

Within the warranty period

Free repairs will be provided during the 1-year warranty period according to what is stated on the Warranty Certificate. Provide MEGLY Customer Support with this product, its Warranty Certificate, and documents with the date of purchase clearly stated (such as the invoice) to order repair service.

After the 1-year warranty has ended

For repairs after the end of the warranty period, consult MEGLY Customer Support. Repairs will be provided for a fee.

1. If this product fails within the warranty period when used correctly in accordance with the precautions in the Instruction Manual, repairs will be provided free of charge (excluding accessories and consumables).
 - a) When requesting a free repair, provide MEGLY Customer Support with this product, its Warranty Certificate, and documents with the date of purchase clearly stated (such as the invoice) to order free repair service.
 2. Even within the warranty period, repairs will be charged a fee in the following cases:
 - a) Product damage due to incorrect use, inappropriate repair, disassembly, or modification
 - b) Product damage due to transport or dropping after purchase
 - c) Product damage due to fire, earthquake, water damage, lightning, other natural disasters, or pollution
 - d) If the Warranty Certificate and the document with the date of purchase clearly stated (such as the invoice) are not provided
 - e) If the customer name has not been written on the Warranty Certificate, or it has been altered
 - f) Product failure or damage due to use contrary to that indicated in the Instruction Manual
 - g) Product failure or damage due to neglect of cleaning after use
 - h) Product failure or damage that occurs when used other than for ordinary home use (for example, for commercial use)
3. The Warranty Certificate is only valid in the country or area of purchase.
 4. Keep the Warranty Certificate in a safe place as it will not be reissued.
 - * Please understand that the personal information (name, address, and telephone number) that you as the customer enter on the Warranty Certificate form may be used for after-sales services during the warranty period and for safety inspection services afterwards.
 - * The Warranty Certificate promises free repairs during the period and under the conditions stated on it. Therefore, the warranty certificate does not limit the legal rights of the customer with regard to the party issuing the warranty certificate (the party responsible for the warranty certificate) or other businesses.

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